

Determining Transgene Copy Number Using Real-Time qPCR on the MJ Research® Opticon™ 2 Continuous Fluorescence Detection System

George Rudenko, Ph.D.¹, Richard Kurtz, Ph.D.², David Batey, Ph.D.², and Virginia Walbot Ph.D.¹

¹Dept. of Biological Sciences, Stanford Univ. ²MJ Research, So. San Francisco, CA

Real-Time Detection

Application Note Vol.2, No.11

Abstract

Real-time quantitative PCR (qPCR) was used to investigate differences in transgene copy number from multiple transgenic lines of maize (*Zea mays L.*). Detection of amplified template was accomplished with SYBR Green I (SGI) chemistry using an MJ Research DNA Engine Opticon® 2 fluorescence detection system. Transgene levels varying from 2 to 21 copies per genome were detected among the lines, demonstrating the sensitivity and target-specificity of the Opticon 2 system for template quantification.

Introduction

The *MuDR/Mu* (*Mu*) transposable elements comprise a complex family of maize transposons that were originally identified in a *Mutator* line of maize with an exceptionally high mutation rate.¹ *Mutator* activity is controlled by the *MuDR* transposon, which consists of the *mudrA* transposase and *mudrB* helper genes, both of which are predicted to encode several distinct proteins as a result of alternative splicing.^{2,3,4,5,6}

The *MuDR/Mu* transposons are the only eukaryotic elements that exhibit two types of transposition. In somatic tissues, *MuDR/Mu* elements exhibit “cut and paste” transposition, in which the element is excised from one site and inserted at a new site. This transposition results in an empty site at the original *Mu* element location. In contrast, just before meiosis and in haploid cells, *MuDR/Mu* elements exhibit “net replicative” transposition, in which the element is copied and inserted into a new site without being excised from the original site.^{7,8}

Transgenic modification of maize is one approach being used to determine which *mudrA* and *mudrB* encoded proteins (termed *MURAs* and *MURBs*, respectively) contribute to specific transposition reactions. Maize lines expressing various combinations of these genes are generated by first crossing transgenic plants expressing individual *MURAs* to lines with or without *MURBs*, and then crossing resulting lines to lines expressing reporter alleles. The crossing schemes produce families in which as many as three transgenes must segregate appropriately. Consequently, each plant must be genotyped to verify that it has a *mudrA*, a *mudrB*, and a reporter allele.

In this study, a real-time qPCR protocol has been developed to genotype transgenic plants by determining *mudrA* transgene copy number. Total genomic DNA was used as the template for reactions utilizing DyNAzyme™ II DNA polymerase* in conjunction with the indicating fluorophore SYBR Green I. The accumulation of amplified nucleic acids was monitored by real-time fluorescence detection on the Opticon 2 system.

Materials and Methods

Quantitative PCR

A 20ng sample of DNA from each maize plant was used as the template for qPCR. DyNAzyme II DNA polymerase was obtained from Finnzymes (F-503L), and SGI was from Molecular Probes (S-7567). Reaction components were assembled in white Hard-Shell® microplates with white wells (MJ Research #HSP-9655) and sealed with ultraclear strip caps (MJ Research TCS-0803). Volumes of individual components and final reaction concentrations are listed in Table 1.

Table 1. Quantitative PCR Reaction Components

Component	Stock concentration	Volume	Final Concentration
DNA sample	10mg/ml	2.0µl	1mg/ml
Forward Primer	(10µM)	0.4µl	0.2µM
Reverse Primer	(10µM)	0.4µl	0.2µM
Mg-free buffer	(10x)	2.0µl	1x
MgCl ₂	(50mM)	0.8µl	2mM
dNTPs	(6.25mM each)	0.4µl	125µM each
SYBR Green I	(10x)	1.0µl	.5x
DyNAzyme II	(2U/µl)	0.2µl	20U/ml
ddH ₂ O		12.8µl	
Total Volume		20 µl	

PCR primers were designed such that the forward primer hybridized within the S-epitope tag nucleotide sequence at the 5' end of the reading frame and the reverse primer hybridized within the *mudrA* coding sequence (see Figure 1). The primer set generated a 308bp amplicon: forward, 5'-GAG-ACC-GCC-GCG-GCC-AAG-T-3', and reverse, 5'-TCA-TCC-AAT-ACA-TCT-GAT-ACT-AT-3'.

Following reaction assembly, plates were transferred to the DNA Engine Opticon 2 real-time detection system (MJ Research). Amplification was performed according to the program listed in Table 2. The primer annealing temperature was experimentally determined using the temperature gradient feature of the DNA Engine Opticon 2 system (data not shown). The temperature at which fluorescence was

*This product is sold under licensing agreements with F. Hoffmann-LaRoche Ltd., Roche Molecular Systems, Inc. and the Applied Biosystems Group of Applied Biosystems Corporation. The purchase of this product is accompanied by a limited license to use it in the Polymerase Chain Reaction (PCR) process in conjunction with a thermal cycler whose use in the automated performance of the PCR process is covered by the up-front licensing fee, either by payment to Applied Biosystems or as purchased, i.e., an authorized thermal cycler.

Real-Time Detection

Table 2. Quantitative PCR Program for Transgene Template Amplification

1. 94°C, 3min
 2. 94°C, 15sec
 3. 63°C, 20sec
 4. 72°C, 30sec
 5. 85°C, 1sec
 6. Plate read
 7. Go to step 2, 42 more times
 8. 72°C, 10min
 9. Melting curve analysis: 65°C to 98°C, 0.2°C/read, 1sec hold
 10. 72°C, 10min
 11. 10°C, Forever
- END

measured during each cycle was set below the T_m of the specific amplicon.¹¹

A dilution series of a plasmid harboring the transgene was mixed with 20ng maize genomic DNA that lacked the transgene. This method controlled for PCR-inhibiting compounds that were residual in the maize genomic DNA preparations. Analysis of these samples was used to generate a standard curve of the log of template quantity versus C(t) cycle.

Table 3. Mean and standard deviation of C(t) values for three replicate measurements of one sample from each of ten maize transgenic lines.

Sample	Mean C(t) ± Standard Deviation
1.	31.22 ± 0.12
2.	33.10 ± 0.30
3.	33.51 ± 0.26
4.	31.60 ± 0.35
5.	33.72 ± 0.41
6.	31.28 ± 0.20
7.	35.60 ± 0.36
8.	31.85 ± 0.26
9.	34.63 ± 0.43
10.	30.69 ± 0.19

The real-time assay was precise and reproducible, as indicated by the low variation in C(t) values across replicates of the same sample template. Table 3 lists the mean ± standard deviation of C(t) values for three replicates of one sample from each transgenic line. The standard deviation ranged from 0.12 to 0.43.

The standard curve generated from duplicate samples of a dilution series of the plasmid harboring the transgene is presented in Figure 2. The concentration of the plasmid ranged from 1,200 initial copies to 1.2×10^6 initial copies. For each reaction, plasmid DNA was mixed with 730 genomic equivalents of maize DNA.

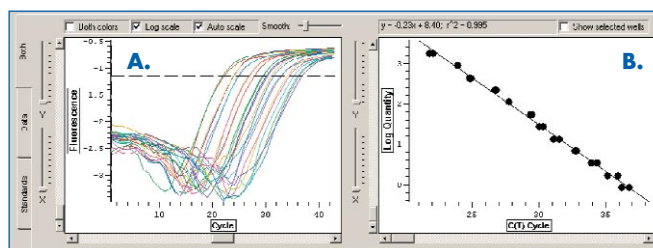


Figure 2. Standard curve for the transgene amplicon, generated from a plasmid serial dilution. A. Quantitation curves for the dilution series, which ranged from 1,200 copies to 1.2×10^6 copies of initial template per reaction. Duplicate reactions were performed at each starting concentration. B. Standard curve plot of log copy number vs. C(t) value. Square of correlation coefficient, $r^2 = 0.995$.

Average transgene copy number per genomic equivalent for each maize line are presented in Table 4. The transgene was present at various levels in the different samples, with a range of approximately 10-fold, from 2 to 21 copies/genome.

Table 4. Average number of transgene insertions in each maize transgenic line.

Sample	Transgene copies/genome (mean ± STD)
1.	16 ± 1
2.	6 ± 1
3.	5 ± 1
4.	13 ± 2
5.	5 ± 1
6.	15 ± 1
7.	2 ± 1
8.	11 ± 2
9.	3 ± 1
10.	21 ± 2

Data Collection and Analysis

The cycle threshold C(t) line was set using the signal/noise ratio option in the Opticon Monitor™ software (MJ Research). The setting used was 10 standard deviations above the mean fluorescence value of the first 3–7 cycles for each sample or standard. This threshold was automatically applied, which ensured consistent analysis of individual samples.

Opticon Monitor software automatically constructed the linear standard curve and determined the transgene copy numbers in the maize samples by interpolation. The maize genome size was estimated to be 2500Mb, and template equivalents of plasmid DNA standards were calculated relative to this genome complexity. All transgenic plants in this study arose from the same parental line and, as such, are expected to have equivalent genome complexities. All samples were analyzed in triplicate within the same plate.

Results

Quantification of Transgene Amounts Using Real-Time qPCR

Real-time qPCR with SGI was performed to investigate *mudrA* transgene copy number in 10 transgenic lines containing a *mudrA* construct. All of these lines were transformed with the same plasmid DNA preparation by the Plant Transformation Facility at Iowa State University (<http://www.agron.iastate.edu/ptf/Web/mainframe.html>). A schematic of the transgene construct is shown in Figure 1.

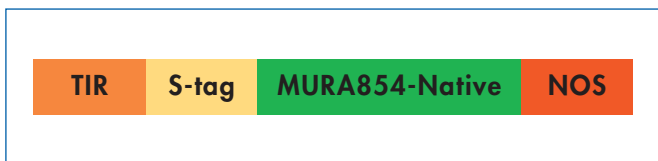


Figure 1. Schematic representation of the *mudrA* transgene (not to scale). The transgene is transcribed from the native promoter (TIR), and is engineered to incorporate a unique S-protein tag (S-tag) to serve in various downstream applications. The transcription terminator (NOS) is derived from the nopaline synthase gene.

Melting Curve Analysis

To assess reaction specificity and to verify product identity, melting curve analysis was performed following amplification. Figure 3 shows the melting curves for triplicate samples of the transgene from one representative maize line. The negative first derivative calculation shows a single peak at approximately 91.5°C, which represents the melting point (T_m) of the amplicon. The presence of a single peak at the same T_m seen in the standards (standard melting curve not shown) supports specific amplification of the transgene.

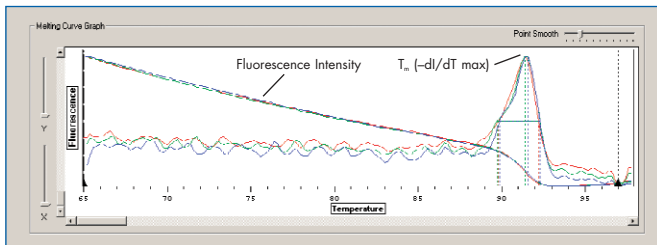


Figure 3. Melting curves depicting the fluorescence intensity and negative first derivative plots for the amplified transgene construct from triplicate measurements of a single sample from a maize plant line.

Discussion

In this study, we developed a real-time qPCR assay to determine the copy number of the *mudrA* transgene in several transgenic maize lines. Results showed that the transgene level varied from 2 to 21 copies per genome among the 10 independent lines. This level of variation is expected, because the use of particle bombardment for transformation of cultured maize embryos is likely to produce multi-copy arrays. In fact, complex transgene arrays containing multiple copies of the transgene plasmid have been previously reported after maize particle bombardment.⁵

Advantages of real-time qPCR for Transgene Quantification

Real-time qPCR can be used to detect and measure the absolute copy number of a transgene in the original lines derived for each transgene construct. Currently, a major limitation in maize experiments is finding individual carriers after completing complex crossing schemes. Such schemes can involve hundreds of plants and varying copy numbers of multiple types of transgene constructs. In some experiments, two discrete loci containing one or more transgenes might be produced. In such instances, the protocol we have developed would permit quantification of the copy number in lines inheriting one or both loci.

The specificity and high signal-to-noise ratio characteristic of real-time qPCR provides several advantages for quantification of transgene copies compared to a technique such as the Southern blot.^{9,10} Compared to real-time qPCR, the DNA blot hybridization technique is costly in terms of reagents, labor, and time: each lane requires 2µg of DNA for radioactive detection or 10µg for fluorescence detection. Consequently, a large tissue sample is required, and so experiments must be postponed until plants have regenerated enough tissue for surveying transgene copy number. Additional complications in determining copy number may exist if rearranged copies are present.

The unique forward primer based on the S-epitope tag sequence and the reverse primer complementary to the *mudrA* gene used in our experiments were designed to hybridize specifically within the transgene coding sequence. This approach should provide superior detection sensitivity by helping to avoid the potential complication caused by the two primers hybridizing to homologous *mudrA* targets present in the maize genome.¹²

The first step in characterizing transformed plants is verifying that the transgene has stably integrated into the chromosome. The second task is to estimate how many copies of the transgene are present and are contributing to the level of transgene expression. Given the relatively long life cycle of maize and the complexity of the genetic analysis, the high degree of accuracy and robustness of the real-time qPCR technique combined with the advantages of high-throughput and reduced assay costs, make this the preferred method over traditional blotting techniques.

Research supported in part by a grant to Virginia Walbot from the National Institutes of Health.

References

1. Robertson, D.S. Characterization of a mutator system in maize. *Mutat Res* 51:21-28 (1978).
2. Hershberger, R.J., Warren, C.A. and Walbot, V. Mutator activity in maize correlates with the presence and expression of the Mu transposable element Mu9. *Proc Natl Acad Sci* 88:10198-10202 (1991).
3. Hershberger, R.J., Benito, M.-I., Hardeman, K.J., Warren, C., Chandler, V.L. and Walbot, V. Characterization of the major transcripts encoded by the regulatory MuDR transposable element of maize. *Genetics* 140:1087-1098 (1995).
4. Eisen, J.A., Benito, M.-I. and Walbot, V. Sequence similarity of putative transposases links the maize Mutator autonomous element and a group of bacterial insertion sequences. *Nucleic Acids Res* 22:2634-2636 (1994).
5. Raizada, M.N. and Walbot, V. The late developmental pattern of Mu transposon excision is conferred by a cauliflower mosaic virus 35S-driven MURA cDNA in transgenic maize. *Plant Cell* 12:5-21 (2000).
6. Lisch, D., Girard, L., Donlin, M. and Freeling, M. Functional analysis of deletion derivatives of the maize transposon MuDR delineates roles for the MURA and MURB proteins. *Genetics* 151:331-341 (1999).
7. Levy, A.A. and Walbot, V. Regulation of the timing of transposable element excision during maize development. *Science* 248:1534-1537 (1990).
8. Raizada, M., Benito, M.-I., and Walbot, V. The MuDR transposon terminal inverted repeat contains a complex plant promoter directing distinct somatic and germinal programs. *Plant J* 25(1):79-91 (2001).
9. Ingham, D.J., Beer, S., Money, S. and Hansen, G. Quantitative real-time PCR assay for determining transgene copy number in transformed plants. *Biotechniques* 31:132-140 (2001).
10. Mason, G., Provero, P., Vaira, A.M. and Accotto, G.P. Estimating the number of integrations in transformed plants by quantitative real-time PCR. *BMC Biotechnology* 2:20 (2002).
11. Optimal plate read temperature for real-time PCR with SYBR Green I. MJ Research Technical Note #104.
12. Rudenko, G. and Walbot, V. Expression and post-transcriptional regulation of maize transposable element MuDR and its derivatives. *Plant Cell* 13:553-570 (2001).

DyNAzyme II Licensing Information

A license under U.S. Patents 4,683,202, 4,683,195 and 4,965,188 or their foreign counterparts, owned by Hoffmann-La Roche Inc. and F. Hoffmann-La Roche Ltd ("Roche"), has an up-front fee component and a running-royalty component. The purchase price of this product includes limited, non-transferable rights under the running-royalty component to use only this amount of the product to practice the Polymerase Chain Reaction ("PCR") and the related processes described in said patents solely for the research and development activities of the purchaser when this product is used in conjunction with a thermal cycler whose use is covered by the up-front fee component. Rights to the up-front fee component must be obtained by the end user in order to have a complete license. These rights under the up-front fee component may be purchased from Perkin-Elmer or obtained by purchasing an Authorized Thermal Cycler. No rights to perform or offer commercial services of any kind using PCR, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is hereby granted by implication or estoppel. Further information on purchasing licenses to practice the PCR Process may be obtained by contacting the Director of Licensing at the Perkin-Elmer Corporation, 850 Lincoln Centre Drive, Foster City, California 94404 or at Roche Molecular Systems, Inc., 1145 Atlantic Avenue, Alameda, California 94501.

Practice of the patented polymerase chain reaction (PCR) process requires a license. The DNA Engine Opticon 2 system includes an Authorized Thermal Cycler and may be used with PCR licenses available from Applied Biosystems. Its use with Authorized Reagents also provides a limited PCR license in accordance with the label rights accompanying such reagents. Some applications may also require licenses from other third parties.

DNA Engine Opticon, Hard-Shell, MJ Research and the helix logo, Opticon, and Opticon Monitor are trademarks belonging to Bio-Rad Laboratories, Inc.

DyNAzyme is a trademark of Finnzymes Oy. SYBR is a trademark of Molecular Probes, Incorporated.

MJ Research, Inc. is the authorized distributor of Finnzymes products in the U.S., Canada, and Japan.

Applications Laboratory
MJ Research, Incorporated
7000 Shoreline Court
South San Francisco, CA 94080



(888) 735-8437 • info@mjr.com • www.mjr.com